

# Meditouch

PRODUCT INFORMATION

# Meditouch

PRODUCT INFORMATION

Total wound care management



**ILDONG**

Head office  
2, Baumoe-ro 27-gil, Seocho-gu, Seoul, Korea Tel : +82-2-526-3114 Fax : +82-2-526-3031  
Factory  
Finished products : 25, Gongdan 1-ro, Anseong-si, Gyeonggi-do, Korea  
Raw materials : 53, Gongdan-ro 98 beon-gil, Heungdeok-gu, Cheongju-si, Chungcheongbuk-do, Korea  
E-mail address : export@ildong.com Home page : www.ildong.com/eng

**ILDONG** PHARMACEUTICAL CO., LTD.

# Classification by purpose

# Total wound care management


## Burns



## Bedsores



## Minor wounds / Cut

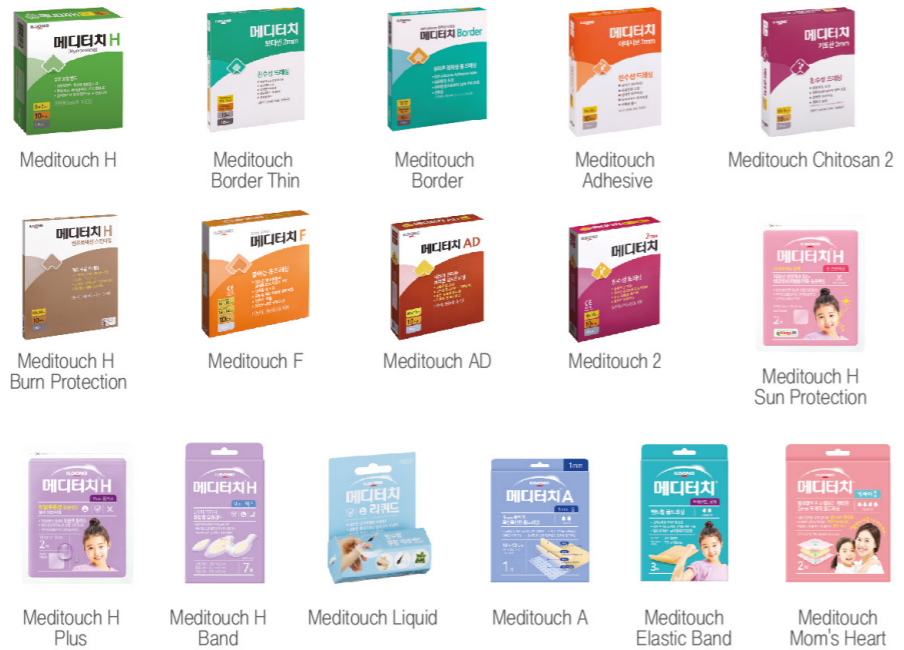

## Surgical wounds



## Foot ulcer / Heel pain



## Finger replantation



## Breast surgery (Procedure)



## Acne / Laser resurfacing



# Meditouch Chitosan 2, 5mm

Meditouch chitosan with high permeability for wounds with excessive exudation and bleeding




## Product features


- Non-adhesive polyurethane foam dressing, containing chitosan
- Meditouch Chitosan 2 (2mm thickness) – Available in 6 sizes
- Meditouch Chitosan 5 (5mm thickness) – Available in 6 sizes

## Point

1. Chitosan-containing foam dressing differentiated from the existing foam dressing
2. Released in large size to reflect market needs
  - 2mm thickness: 10x20cm, 20x20cm
  - 5mm thickness: 35x50cm, 30x100cm, 35x150cm
3. Prevention of bedsores during surgery, management of bedsores in ICU wards, applied to wounds with excessive exudation and bleeding
4. Cut and use according to the size of the wound (zipline pouch applied, except for small sized products)

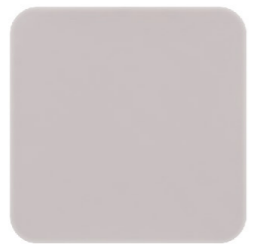


Meditouch Chitosan 2, 5mm  
(10x10 cm)

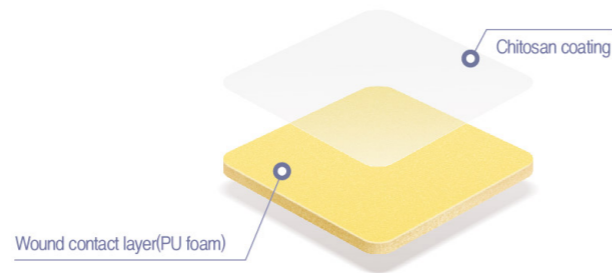

Product	Size (cm)	Pieces/Box	Product	Size (cm)	Pieces/Box
Meditouch Chitosan 2	5 x 10	10	Meditouch Chitosan 5	10 x 10	10
	5 x 20	10		10 x 20	10
	5 x 30	10		20 x 20	10

# Meditouch Chitosan Roll 7mm

Meditouch chitosan with high permeability for wounds with excessive exudation and bleeding

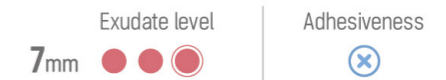


## Product features

- Non-adhesive polyurethane foam dressing containing chitosan
- 7mm thickness, the thickest roll type
- A roll type, easy to carry compared to other large-sized types
- Available in 2 sizes

## Point

1. Chitosan-containing foam dressing differentiated from the existing foam dressing
2. Released in large size to reflect market needs – 20x100cm
3. Increased absorbency with 7mm thickness foam compared to competing products (thickness of competing products are 5mm)
4. For exudate, excessive bleeding wounds, wounds that are difficult to apply with the existing foam sizes, management of operating wounds and bedsores
5. Cut and use according to the size of the wound (zipline pouch applied)

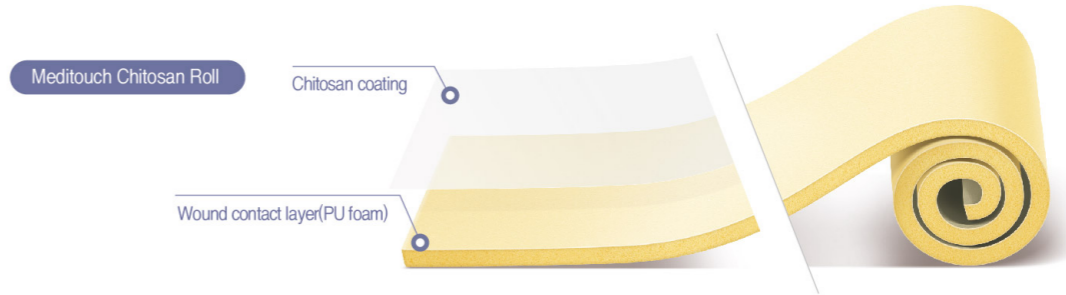

Product	Size (cm)	Pieces/Box
Meditouch Chitosan Roll	10 x 65	1
	20 x 100	1

# Meditouch Heel S/L

8mm

Non-adhesive polyurethane foam dressing for feet  
Suitable for bedsores, diabetic feet, foot ulcers and splints

Exudate level ●●● Adhesiveness ⊗



## Product features

- Non-adhesive polyurethane foam dressing for feet
- **The thickest 8mm foam dressing** compared to competing products
- Released in 2 sizes (small and large size) **which can be selected and applied according to the foot area**
  - ▶ Small – covers heel
  - ▶ Large – covers heel and around the ankle

## Point

1. For prevention and management of **bedsores on the heel and ankle of ICU immobilized patients**
2. Used as a cushion for **heel pain relief** and bed sore prevention **after OS foot cast**
3. **The thickest 8mm foam dressing** compared to competing products **absorbs exudate and provides excellent cushioning**

Meditouch Heel S



Meditouch Heel L



Product	Size (cm)	Pieces/Box
Meditouch Heel S	13 x 10	5
Meditouch Heel L	22 x 11.5	5

# Meditouch Net

5mm

A **perforated net layer**(PU film) to the 5mm thickness of absorbing layer helps to **maintain its shape** when absorbing excess exudates

Exudate level ●●○ Adhesiveness ⊗



## Product features


- A **perforated net film** under the absorbing layer helps to maintain a shape of foam dressing
- Uses for wounds with **heavy exudates**
- **No residue left** around the wound when changing the dressing

## Point

1. A **perforated layer** can absorb wound exudates, results in **longer wear time and reducing risk of maceration**
2. A perforated layer allowing free passage of wound exudates to the foam **prevents skin maceration**
3. Easier and convenient use for patients providing self-care at home




Meditouch Net 5mm (10x10cm)

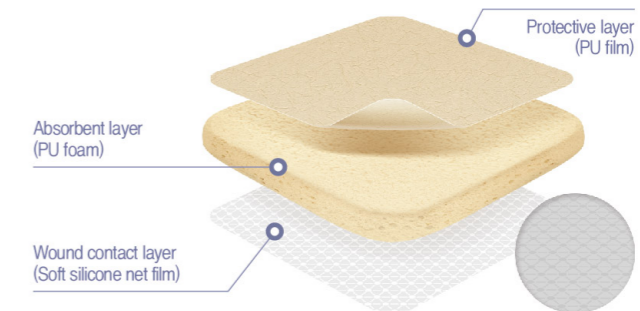

Product	Size (cm)	Pieces/Box
Meditouch Net 5mm	10 x 10	10
	10 x 20	10
	20 x 20	10
	30 x 55	10

# Meditouch 2, 5mm

Highly absorbent polyurethane foam dressing with 3 layers and thickness of 2mm and 5mm

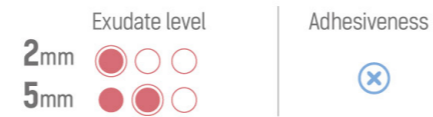

## Product features


- A polyurethane foam dressing with 3 layers and thickness of 2mm and 5mm
- A thickness of 2mm to little exudates and 5mm to excessive exudates
- A secondary dressing required

## Point

1. An absorbent polyurethane foam dressing supporting moist wound treatment
2. A high density of foam prevents leakage of wound exudates when the absorption layer is placed under pressure



Meditouch 2, 5mm  
(10x10cm)

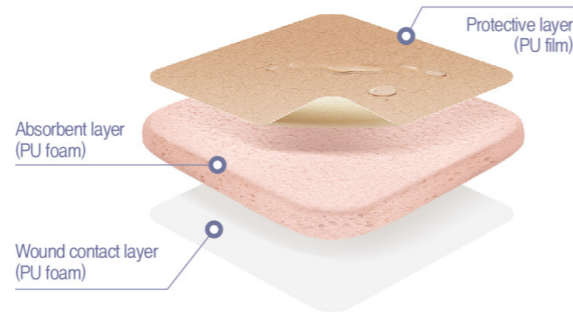

Product	Size (cm)	Pieces/Box	Product	Size (cm)	Pieces/Box	Size (cm)	Pieces/Box
Meditouch 2	5 x 10	10	Meditouch 5	10 x 10	10	35 x 50	10
	10 x 10	10		10 x 20	10	30 x 100	5
	5 x 20	10		20 x 20	10	35 x 150	5

# Meditouch Hand

The only wet dressing in Korea exclusively designed for hands wound care



## Product features


- The only wet dressing in Korea exclusively designed for hands wound care
- Contained hand-shaped molding to secure the wound hands increases the convenience of treatment/use for healthcare professionals
- Available in 4 sizes

## Point

1. Designed for injured fingers and lower part of the hands to contact easily to the absorbent layer
2. Ease of use for healthcare professionals and patients



One touch hand dressing

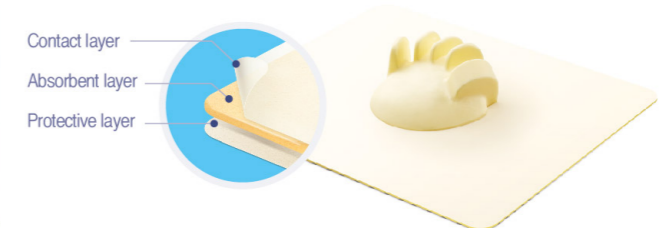

Product	Size (cm)	Pieces/Box
Meditouch Hand S	Hand 9 x 7.5 x 5 Pad 31 x 36	1
Meditouch Hand L	Hand 12.5 x 13 x 7.5 Pad 31 x 36	1

# Meditouch Border Thin

2mm

A convenient and light product, Meditouch Border Thin can be used in exposed body areas like face, neck and arm

Exudate level ● ○ ○ | Adhesiveness ●



## Product features


- Soft silicone border type adhesive dressing that does not require secondary dressing
- 2mm thickness foam
- No peeling paper
- Ivory PU film
- Available in 6 sizes

## Point

1. Easy to apply to various wounds with 6 sizes
2. Appropriate for the management of exposed wounds due to its thin thickness
3. Passed European standard for waterproof test



Meditouch Border Thin  
5.5x6.5cm (foam 3x4cm)

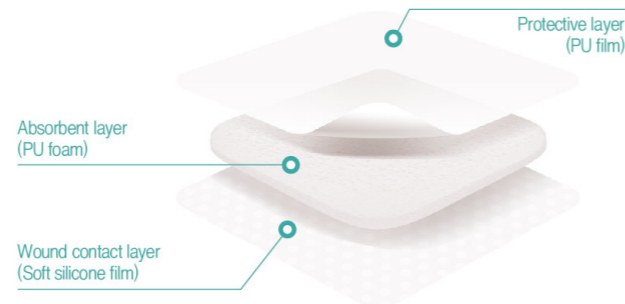

Size (cm)	Pieces/Box
5.5 x 6.5 (foam 3 x 4)	30
9 x 9 (foam 5 x 5)	10
10 x 10 (foam 7 x 7)	10
9 x 15 (foam 5 x 10)	10
9 x 25 (foam 5 x 20)	10
15 x 15 (foam 10 x 10)	10

# Meditouch Adhesive

2mm

Increased adhesion, minimized irritation  
Soft silicone adhesive type with large size, Meditouch Adhesive!

Exudate level ● ○ ○ | Adhesiveness ●



## Product features


- 2mm thickness of adhesive polyurethane foam dressing
- Soft silicone adhesive dressing
- Silicone gel is coated on polyurethane layer
- Available in 4 sizes

## Point

1. Soft silicone adhesive method increases adhesion and minimizes irritation
2. Does not stick to the wound surface when removing dressing
3. Released in 35x50cm large size (Reflecting market needs/Burn)



Meditouch Adhesive  
(10x10cm)

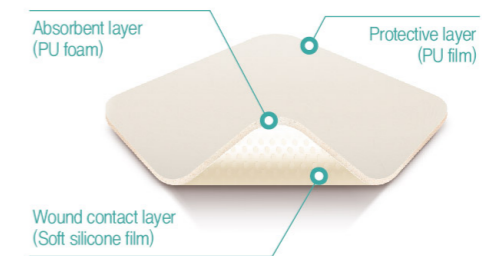

Size (cm)	Pieces/Box
10 x 10	10
10 x 20	10
20 x 20	10
35 x 50	10

# Meditouch AD

2mm

A polyurethane foam dressing that has an adhesive **soft silicone layer at the wound surface**. An adhesiveness helps patients to practice wound self-care at home.




## Product features

- An adhesive polyurethane foam dressing with a thickness of 2mm
- Available in 6 sizes
- A **soft silicone layer at the wound contact surface** makes it to attach easier
- Meditouch AD Circle has a thickness of 12 mm and **special round shape** to apply after breast surgery

## Point


1. A **self-adhesiveness** on the wound contact layer **provides convenience** to use for healthcare professionals and patients
2. Although a secondary dressing is required, a self-adhesiveness at the wound contact surface helps patients to practice **self-dressing changes at home easily**
3. Meditouch AD Circle is specially designed **for breast** and **to fix drain tubes**

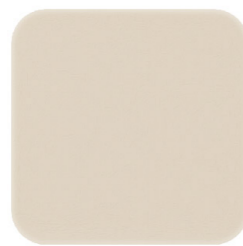

Meditouch AD (10x10cm)



Meditouch AD Circle

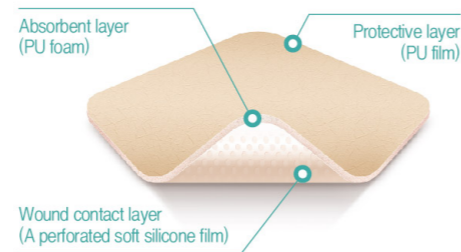

Size (cm)	Pieces/Box	Size (cm)	Pieces/Box
5 x 10	10	10 x 20	10
5 x 20	10	20 x 20	10
10 x 10	10	12mm Circle 6.5 x 1.6	20

# Meditouch H & Meditouch H Sun Protection

A hydrocolloid dressing for light exudates, Used for light wounds on the face, acne, and after laser treatments



## Product features


### Meditouch H Sun protection

- Acquired **Allergy Friendly Product Award certificate**
- US FDA registered
- Sheet type: Applied 'Skin Wrinkle Processing' technology
- Dot type: Applied **beveled edge**
- UPF50+ scar protection from the sun → **Prevents skin pigmentation**

## Point

1. A hydrocolloid dressing **for light wounds**
2. **Better convenience** than foam dressing as **no secondary dressing is required**
3. Available in various sizes for healthcare professionals and patients
4. (Meditouch H Sun protection) Beveled edges help the dressing border to stay flat and smooth on the patient's skin, **reduce the risk of pushed off when putting on or taking off clothes, and allow to wear a light makeup**

### Dot type

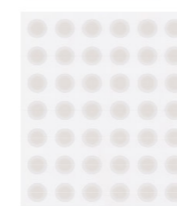


Meditouch H Sun Protection (dot 10mm)

### Sheet type



Meditouch H Sun Protection (Sheet 10x10cm)

### Band type

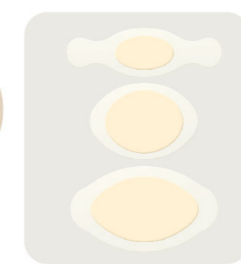

Product	Size (cm)	Pieces/Box	Product	Size (cm)	Pieces/Box	Product	Size (cm)	Pieces/Box
Meditouch H	dot 9mm	42R x 50	Meditouch H Sun Protection	dot 10mm	42R x 50	Meditouch H Band (Mixed type)	S (2 x 6.3) 3EA M (3.1 x 4.4) 2EA L (3.5 x 6) 2EA	
	10 x 10	10		dot 12mm	42R x 50			
Meditouch H Sun Protection	10 x 10	10	Meditouch H Sun Protection (Mixed type)	Dot 12mm x 30R				
	20 x 20	10		Dot 10mm x 18R				
	10.5 x 10.5	2		Dot 12mm x 30R				
				Dot 10mm x 36R				

# Meditouch Border

**3mm**

An adhesive soft silicone dressing applied to the wound contact layer **does not require a secondary dressing** and **minimizes skin trauma** during dressing removal

Exudate level ● ◐ ◑ | Adhesiveness ●




## Product features

- A polyurethane foam dressing with a thickness of 3mm
- Available in 6 sizes
- An adhesive soft silicone helps to maintain a moist wound environment, and reduces skin trauma and pain during dressing removal
- High absorption and retention

## Point

1. No secondary dressing required as soft silicone film used at the wound contact layer
2. A high tack makes the foam stay over curves and irregular surfaces with ease
3. Can be reattached after checking the wound area

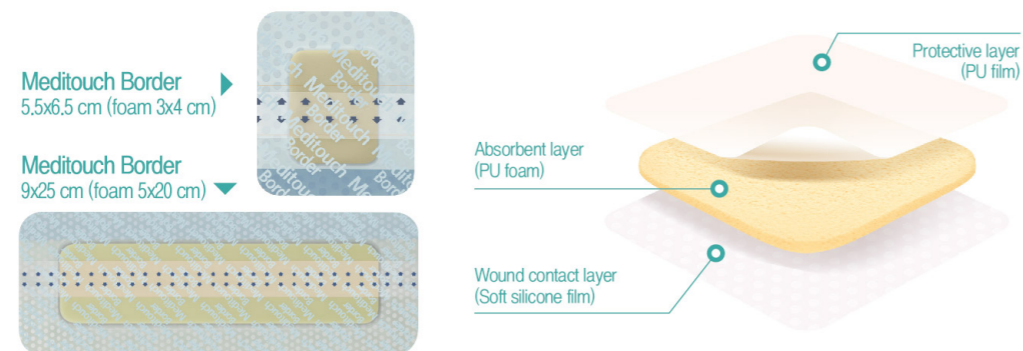

Size (cm)	Pieces/Box	Size (cm)	Pieces/Box
5.5 x 6.5 (foam 3 x 4)	30	9 x 15 (foam 5 x 10)	10
7.5 x 7.5 (foam 5 x 5)	10	9 x 25 (foam 5 x 20)	10
10 x 10 (foam 6.5 x 6.5)	10	15 x 15 (foam 10 x 10)	10

# Meditouch Border Post-op

**1.7mm**

A product with high absorbency after surgery and minimizing skin irritation and pain during dressing change

Exudate level ● ● ● | Adhesiveness ●




## Product features

- A soft silicone contact layer enables to remove easily and minimizes skin irritation and pain when dressing removal
- Made of high absorbent non-woven fabric absorbs exudates approximately 1.8 times more than polyurethane foam dressing, doubles retention capacity of exudates
- Bordered shape fits onto curved and contoured areas and has great adhesion to skin
- Available in 6 sizes

## Point

1. A synthetic gauze dressing for the protection of sutures after surgery
2. Applied 'Flexible Cut' to the absorbent layer allows wound areas move flexibly while wearing, even allowing for light exercise
3. The waterproof top barrier allows patients to wash or shower

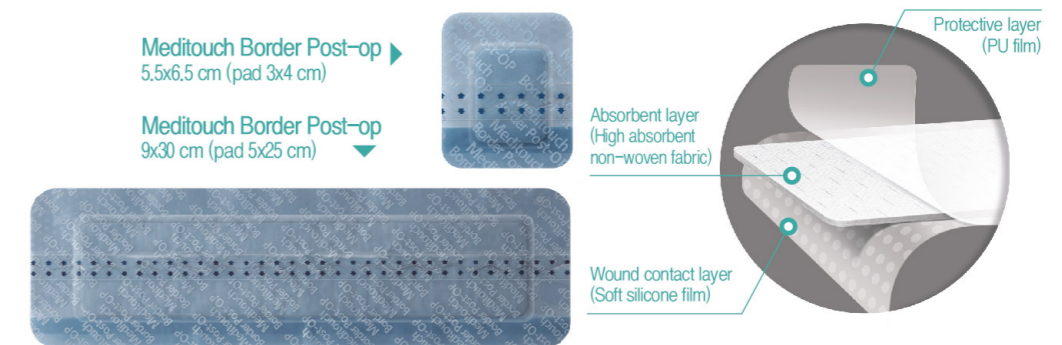

Size (cm)	Pieces/Box	Size (cm)	Pieces/Box
5.5 x 6.5 (pad 3 x 4)	30	9 x 20 (pad 5 x 15)	10
7.5 x 7.5 (pad 5 x 5)	10	9 x 25 (pad 5 x 20)	10
9 x 15 (pad 5 x 10)	10	9 x 30 (pad 5 x 25)	10





# Meditouch F

2mm

A thickness of 2 mm polyurethane foam dressing with adhesive acrylic film on the top, no required secondary dressing




## Product features


- Applied adhesive acrylic film on the surface of polyurethane foam
- Available in 5 sizes
- No required fixation tape or secondary dressing

## Point

1. Applying an adhesive acrylic film on top of the absorbent layer(PU foam)
2. No required a secondary dressing provides higher convenience for healthcare professionals and patients

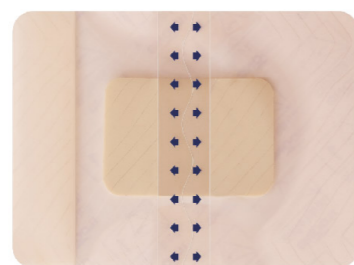

Meditouch F 2mm  
7.5x10 cm (foam 3.5x5 cm)

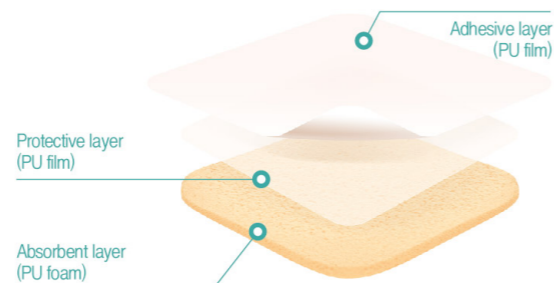

Size (cm)	Pieces/Box
7.5 x 10 (foam 3.5 x 5)	50
9 x 15 (foam 5 x 10)	10
14 x 14 (foam 10 x 10)	10
9 x 25 (foam 5 x 20)	10
9 x 35 (foam 5 x 30)	10



# Meditouch A

1, 5mm

A medical silicone adhesive foam dressing to excess exudates and for patients with sensitive skin



## Product features

- **High exudate absorption:** Contained Sodium Alginate absorbs exudates upto 15 to 20 times more than its initial molecular weight
- **Sodium Alginate** forms a protective film on the wound and prevents bacterial infection at the wound site
- **Hemostatic effect** from ion exchange occurs as alginate absorbs wound exudate
- The medical soft silicone foam to a wound contact layer makes it **hypoallergenic**

## Point

1. For patients to wash or shower, a **waterproof non-woven fabric** included
2. **Available wounds**
  - An open wounds such as **abrasion** from falls, a **laceration** from a deep cut or tearing of your skin by sharp items, and **burns** which has **heavy bleeding and exudates**
  - Wounds **after surgery(stitches)**, **bedsores** with heavy exudate
3. **Recommended** a thickness of 5mm for the first 1 or 2 days after wound and 1mm from the 3rd day



Product	Size (cm)	Pieces/Box
Meditouch A 1mm	4 x 9	3
	10 x 10	1
Meditouch A 5mm	10 x 10	1

# Meditouch Burn Protection

## 2, 5mm

A polyurethane foam dressing, containing a complex of silver and activated carbon, used for wounds with a large amount of exudates such as burns

Exudate level

2mm ●○○○

5mm ●●○○

Adhesiveness

⊗

### Product features

- A thickness of 2mm & 5mm foam dressing (A fixation tape included)
- Silver nitrate contained in a wound contact layer has antibacterial activity
- An activated carbon contained in a wound contact layer reduces wound odor by absorbing gases released by bacteria
- High absorption and retention capacity as PU foam has pores and cells

### Point

1. For patients with minor burns or those who find it difficult to visit clinics/hospitals
2. Recommended to patients looking for the product specially designed for burn wound care

#### TOP PU film backing


Maintaining a moist environment for wound healing and preventing contamination of the polyurethane foam

#### PU foam absorbent layer

A wound contact layer containing a complex of silver and activated carbon, absorbs wound exudate and creates a moist environment

#### OPP film

A release film preventing contamination of the wound contact surface until the product is used

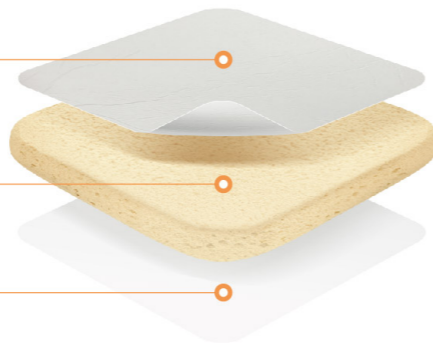

Product	Size (cm)	Pieces/Box
Meditouch Burn Protection 2mm	foam 5 x 7 (fixation tape 10 x 12)	3
Meditouch Burn Protection 5mm	foam 10 x 10 (fixation tape 15 x 15)	1

# Meditouch Elastic Band

A waterproof foam dressing, thin and highly elastic for comfortable movement

Exudate level


●○○○

Adhesiveness

●

### Product features & Point

- No secondary dressing required
- A polyurethane foam dressing providing moist environment for wound healing
- Used waterproof non-woven fabric for patients to wash or shower
- The thin and elastic film can offer comfort during movement

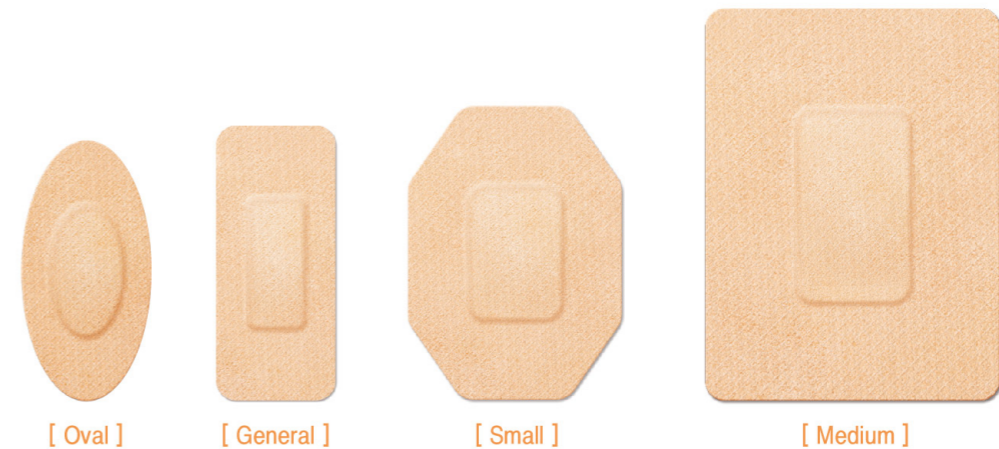

Product	Size (cm)	Pieces/Box
Meditouch Elastic band (S)	5.5 x 7.5 (foam 2.5 x 3.5)	3
Meditouch Elastic band (M)	7.5 x 10 (foam 3 x 5)	3
Meditouch Elastic band (General)	3 x 7 (foam 1.5 x 3.4)	7
Meditouch Elastic band (Mix)	General 3 sheets, S 2 sheets, Oval 3 sheets(foam 2 x 3.5, film 3.5 x 7)	